


Wheelchair lifts

Autoadapt | How to comfortably and safely enter your vehicle seated in a wheelchair.





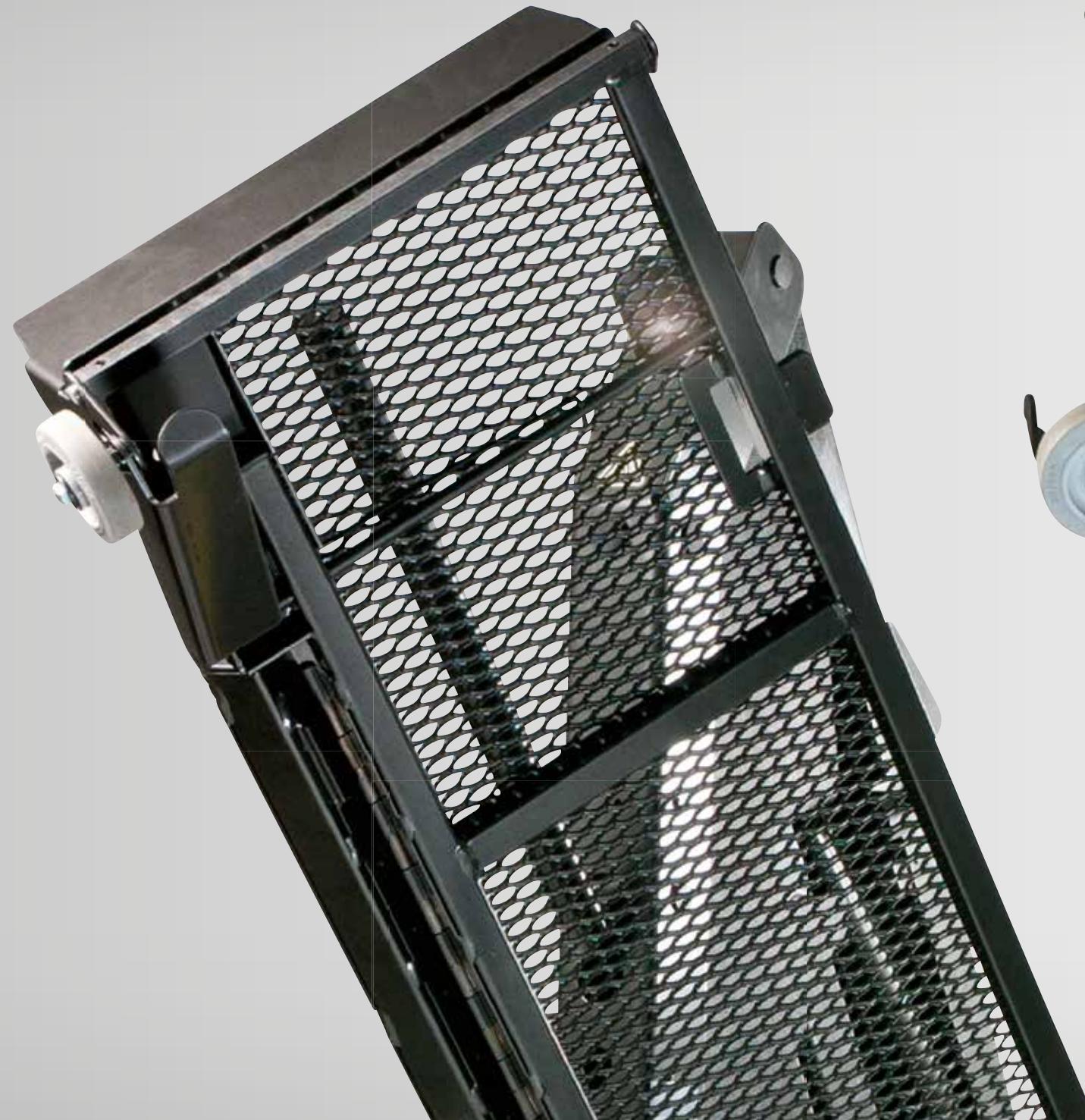

Welcome to
Autoadapt,
we hope you'll find this to be
an uplifting experience

THE SAFEST WAY for an adult person to travel in any car is seated in a car seat wearing a seat belt. We strongly believe that the first rule of vehicle adaptation is to never compromise the level of safety. With the proper protection it is possible to enjoy the ride seated in a wheelchair without compromising safety.

Over 40 years ago Ralph Braun had the need to enter his vehicle seated in his tri-wheeler. He quickly realized he needed some sort of lift and since no solution was available at that time he had to invent one himself. Since then our part-owner BraunAbility has been perfecting their wheelchair lifts and today they are the largest manufacturer of wheelchair lifts in the world, well respected for their quality and foremost their reliability.

The following pages will introduce you to a number of different types of lifts that all will help you and your wheelchair get into the vehicle. However in order to find the most suitable solution for you we encourage you to contact one of our authorised dealers for an assessment.



Platform lighting ensures visibility in low light operation and makes the user feel safer.

Quick facts Millennium

- A strong, solid lifting platform featuring both front and side entry
- Unique roll stop. Powered by hydraulics and locked mechanically, providing additional safety
- Illuminated platform
- Bridging feature simplifies boarding from pavement
- Durable hand-held control box
- Rigid handrail for wheelchair users and standees
- Rigid vertical mechanical In Board barrier for extra safety and bumper clearance
- Designed for side or rear door installation
- Cycle counter
- EMC test approved
- CE-marked

Does your parking situation allow you to be generous with space around your vehicle?



IF YOU HAVE access to wide parking spaces or don't have to worry about someone blocking the back of your vehicle, just about any BraunAbility wheelchair lift is a good choice. We have platform lifts that suit many types of situations for example the extremely stable Millennium or the telescopic Century XT.

“ A car of your own brings freedom. You don't have to be dependent of others, you can visit friends whenever you want to and essentially mind yourself.

– Håkan Thörnqvist, Sweden

Quick facts Century XT

- Extra long telescopic stacking platform accommodates full-sized mobility scooters
- Illuminated platform
- Bridging feature simplifies boarding from pavement
- Durable hand-held control box
- Rigid handrails for wheelchair users and standees
- Rigid vertical mechanical In Board barrier for extra safety and bumper clearance
- Designed for side or rear door installation
- Cycle counter
- EMC test approved
- CE-marked

Is your parking situation spacious ...
... but free passage and line of sight is important to you?



HAVING A FREE line of sight through the window improves the quality of the ride but it could just as well be a safety issue.

Another thing to consider is free passage. Having the ability to “walk through” the lift in its stowed position could be a real timesaver if you have other passengers hopping on and off. The Vista Split does not only provide a clear view but it is also out of the way when not in use.



IF YOU DON'T want your wheelchair lift to block the view or the passage through the doors when you aren't using it consider a Vista series wheelchair lift, the UVL or the Vangater to fulfil your needs.

The UVL (Under vehicle lift) is as the name implies located under your vehicle where it neither takes up valuable floor space nor blocks any views or exits.



When not in use the Vista Split platform is stowed in parallel to the vehicle's sides, leaving an unobstructed view and good access.



The innovative trifold platform Vangater stows in parallel to the vehicle's door leaving plenty of room for both view and access.

◀ Quick facts Vista Split

- A clear path for both visibility and access
- Illuminated platform
- Bridging feature simplifies boarding from pavement
- Durable hand-held control box
- Cycle counter
- EMC test approved
- CE-marked

Quick facts UVL

- Discreet and out of the way when not in use
- Durable handheld control box including door opener button
- Slip-resistant platform surface
- Lift will not raise unless the automatic roll stop is up

Quick facts Vangater

- Trifold platform featuring front and side entry
- Durable hand-held control box
- Vertical mechanical In Board barrier for extra safety
- Designed for side or rear door installation
- EMC test approved
- CE-marked



“ The best thing about my car adaptation is the sense of freedom to come and go as I please. All my daily chores become so much easier to do.

— Katharina Oskarsson, Sweden

Quick facts Swing-A-Way

- Requires a minimal parking space. The deployed platform protrudes no more from the side of the vehicle than the fully opened front door
- Platform handrail provides added security for wheelchair user
- EMC test approved
- Does not obstruct normal use of door opening
- CE-marked

Does your parking situation require you to parallel park or use tight parking spaces?

IN AN URBAN environment space is a luxury many can't afford. With a lift that protrudes no longer from the side of the vehicle than the open passenger door, tight parking situations and parallel parking is made possible. The lift that does all this is the Swing-A-Way.





Additional products to make the best out of your wheelchair lift



AUTOADAPT 6-WAY BASE

THE AUTOADAPT 6-WAY BASE makes it easier to transfer from a wheelchair to a safe, stable and comfortable car seat inside a van. It allows you to move the driver seat backwards, rotate it and adjust it in height to bring the driver seat next to the wheelchair. The Autoadapt 6-Way Base is crash tested and certified with the female seatbelt buckle integrated in the base. A foot rest that follows the lengthwise movement is available as an option.



More benefits at:
autoadapt.com

WHEELCHAIR RESTRAINTS

HAVING A WHEELCHAIR in a vehicle, occupied or not means you must secure it before driving. Our Unwin wheelchair restraints come in a variety of different types from the four point tie downs to the Dock'N'Lock; all with the sole purpose of safely securing your mobility device. But remember if you travel seated in a wheelchair you must wear a real seat belt and have proper head and neck protection.



Autoadapt worldwide mobility solutions

AUTOADAPT WAS FOUNDED in 1996 and has since grown to one of the world's leading manufacturers of car adaptation solutions for people with reduced or limited mobility. Today Autoadapt is represented in over 50 countries worldwide, with more than 300 trained Autoadapt dealers plus an additional network of about 700 sub-dealers. Autoadapt's head office, development and production plant is located in Stenkullen, in the south western part of Sweden. A subsidiary in Birmingham, Autoadapt UK Ltd., operates on the UK market.

Since June 2013 Autoadapt is part-owner of Unwin Safety Systems located in Somerset, UK. Being a pioneer in safety Unwin has become one of the worlds leading suppliers of wheelchair tie-down and occupant restraint systems. Together with our part-owners BraunAbility and Bruno Independent Living Aids, Autoadapt is the leading global provider of safe vehicle adaptation solutions for people with reduced or limited mobility.



TO GUARANTEE YOUR AND all our end users' safety Autoadapt products are designed and tested according to current directives and standards. By working together with users, organisations, the automotive industry and other interested parties we are constantly pushing the boundaries in order for you to be able to choose the car you want, to travel as safe as possible and live a more mobile and independent life.





OUR AREAS OF EXPERTISE

AT AUTOADAPT WE help people do what so many take for granted; drive to work, pop down to the shops, drop in on friends or go on a great driving holiday. It is with great pride we design and manufacture solutions that people from all over the world put their trust in every day. To help you find the solution to suit your needs, simply look at the corresponding area of expertise.

You are currently reading our brochure on Wheelchair lifts. For a quick look at our other areas of expertise, please see below.

GETTING SEATED

THE SAFEST WAY for an adult person to travel in any car is seated in a car seat wearing a seat belt. To facilitate this we invented a number of solutions that enables you enjoy the ride from the comfort of a car seat, both as passenger and driver.

INDEPENDENT DRIVING

USING FOOT PEDALS to drive is a way that somehow became standard almost a century ago. Fortunately we can change that and enable you to take the wheel and go for a spin.

LIFTING AND STOWING

BEFORE YOU GO for a ride you probably want to make sure you have your wheelchair with you. Regardless of its type and size there is a solution to stow your wheelchair in a safe and secure manner.

WHEELCHAIR RESTRAINTS

OUR PARTNER Unwin Safety Systems has pioneered wheelchair restraints and occupant safety solutions for over 50 years. Unwin continue to introduce innovative design and engineering whilst being at the forefront of international safety standards.



”The wheelchair
enables me to walk,
the car enables me
to run.

– Mona Nordell, Independent driver.

*Autoadapt is the proud supplier of independence to drivers and passengers all over the world.
See their stories at: autoadapt.com*

AUTOADAPT AB
Åkerivägen 7
S-443 61 Stenkullen
Sweden
Phone: +46 302 254 00
E-mail: info@autoadapt.se

AUTOADAPT UK LTD
Unit 1, Windsor Industrial Estate
Rupert Street
Aston, Birmingham B7 4PR, UK
Phone: +44 121 33 35 170
E-mail: info@autoadapt.co.uk

Illustrations, descriptions and specifications in this brochure are based on current product information.
Autoadapt reserves the right for alterations without previous notice. © 2013 Autoadapt AB

 www.autoadapt.com
facebook.com/autoadapt
twitter.com/autoadapt